



## SAFETOUCH™ CATH

STERILE, SINGLE-USE CATHETER  
WITH SAFETY DEVICE FOR IV INFUSION



# SAFETOUCH™ CATH

STERILE, SINGLE-USE CATHETER WITH SAFETY DEVICE FOR IV INFUSION

- Irreversible, passive safety mechanism
- Audible click confirms locking of safety device
- Does not require a change in puncture technique
- Radiopaque
- Ultra sharp three-beveled needle to minimize discomfort
- Siliconization of needle minimizes penetration and gliding force
- Color-coded hub to indicate gauge
- Transparent PUR catheter tube for quick visualization of blood in the flashback chamber to facilitate insertion
- Flashback in hub means the needle is correctly in the vein
- Flashback in catheter means the catheter is correctly in the vein
- Not made with latex, DEHP or PVC
- Available in 16-24 G
- IV infusion duration: hours to several days
- For use by healthcare professionals only

CE 0123

Class IIa Medical Devices

Rule 7

MDD 93/42/EEC

UMDNS 18331

## PRODUCT COMPLIANCE

- CE-marked, Class IIa Medical Devices, Rule 7, MDD 93/42/EEC, UMDNS 18331
- Complies with the following norms, directives, and regulations:
  - ISO 594-1:1986
  - ISO 7864:1993
  - ISO 14644-4:2001
  - ISO 594-2:1998
  - ISO 10555-5: 2013
  - ISO 23908:2011
  - ISO 2859-1:1999 /AMD 2011
  - ISO 10993-10:2010
- Labels contain 20 languages

## MANUFACTURING DETAILS

**Legal manufacturer:** Nipro Corporation

**Country of origin:** Thailand

## STERILIZATION AND SHELF LIFE

### Sterilization

EtO (Ethylene oxide)



### Single-use only



### Indicator sticker


Each outer box (or shipping carton) contains a chemical indicator that indicates sterility. This indicator comes in the form of a [blue] sticker that turns [red] when sterilized.

### Shelf life

5 years



## MATERIALS USED

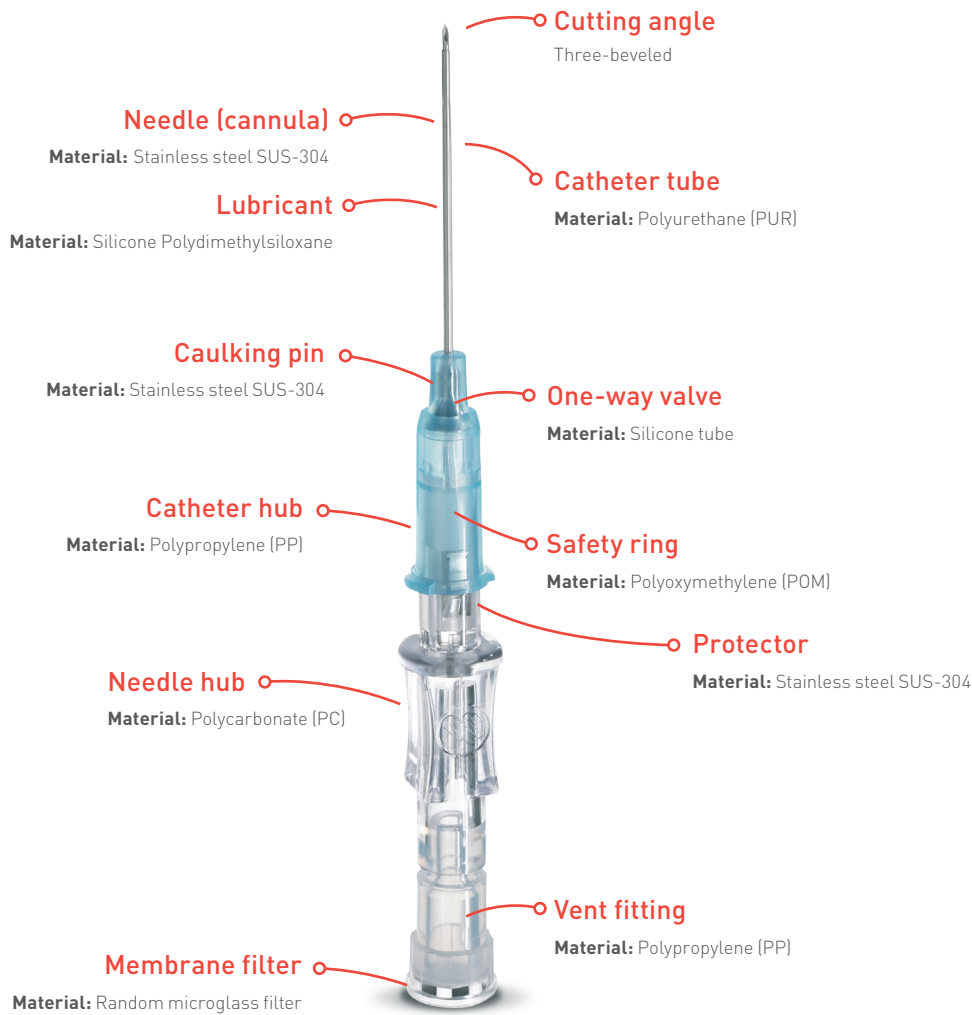

**Not made with  
Latex, DEHP or PVC**

**Glue**  
Material: Acrylic resin

## PRODUCT RANGE OVERVIEW

Reference code	Description	Color	Flow rate not less than
NIC-SF 1650P-ET	IV Cath SF 16Gx2" (1,7x50mm) PUR		≥ 162 ml/min
NIC-SF 1832P-ET	IV Cath SF 18Gx1 1/4 " (1,3x32mm) PUR		≥ 81 ml/min
NIC-SF 1850P-ET	IV Cath SF 18Gx2" (1,3x50mm) PUR		≥ 81 ml/min
NIC-SF 2025P-ET	IV Cath SF 20Gx1" (1,3x50mm) PUR		≥ 49 ml/min
NIC-SF 2032P-ET	IV Cath SF 20Gx1 1/4 " (1,1x32mm) PUR		≥ 49 ml/min
NIC-SF 2050P-ET	IV Cath SF 20Gx2" (1,1x50mm) PUR		≥ 49 ml/min
NIC-SF 2225P-ET	IV Cath SF 22Gx1" (0,8x25mm) PUR		≥ 29 ml/min
NIC-SF 2232P-ET	IV Cath SF 22Gx1 1/4" (0,8x32mm) PU		≥ 29 ml/min
NIC-SF 2419P-ET	IV Cath SF 24Gx3/4" (0,7x19mm) PUR		≥ 17 ml/min
NIC-SF 2425P-ET	IV Cath SF 24Gx1" (0,7x25mm) PUR		≥ 17 ml/min

# PACKAGING DETAILS

Reference code	Gauge	Length	Inner box					Outer box					
	(G)	(mm)	Quantity	L (mm)	W (mm)	H (mm)	Weight (kg)	Quantity	L (cm)	W (cm)	H (cm)	Weight (kg)	m³
NIC-SF-xxxxP-ET	xx	xx	50	160	63	145	9,840	2,000	662	337	319	10,840	0,0712

## Languages

20 languages on inner and outer box:

English (EN), Bulgarian (BG), Czech (CS), Danish (DA), German (DE), Greek (EL), Spanish (ES), Finnish (FI), French (FR), Hungarian (HU), Italian (IT), Lithuanian (LT), Dutch (NL), Norwegian (NO), Polish (PL), Portuguese (PT), Romanian (RO), Russian (RU), Swedish (SV), Turkish (TR)


## Transport conditions







Closed and dry

## Storage conditions

Open the packaging only immediately before use to guarantee sterility.

## Label details

	Individual unit	Inner box	Outer box
Product name	✓	✓	✓
Product size	✓	—	✓
Quantity (pieces)	—	✓	✓
REF	✓	—	—
LOT	✓	—	—
	✓	—	—
STERILE EO	✓	✓	✓
CE	✓	✓	✓

	Individual unit	Inner box	Outer box
	✓	✓	✓
	✓	✓	✓
	✓	✓	✓
	✓	✓	✓
EC REP	✓	✓	✓
	—	✓	✓
	✓	✓	✓

