



PRODUCT DATA SHEET

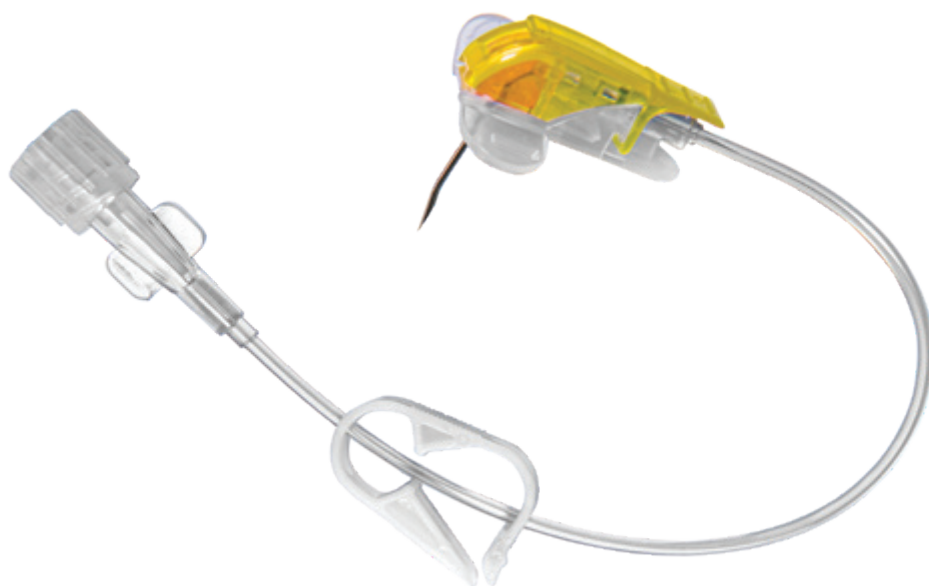

SAFETOUCH HUBER NEEDLE SET

STERILE, SINGLE-USE, NON-CORING HUBER NEEDLE
INFUSION SET WITH SAFETY DEVICE

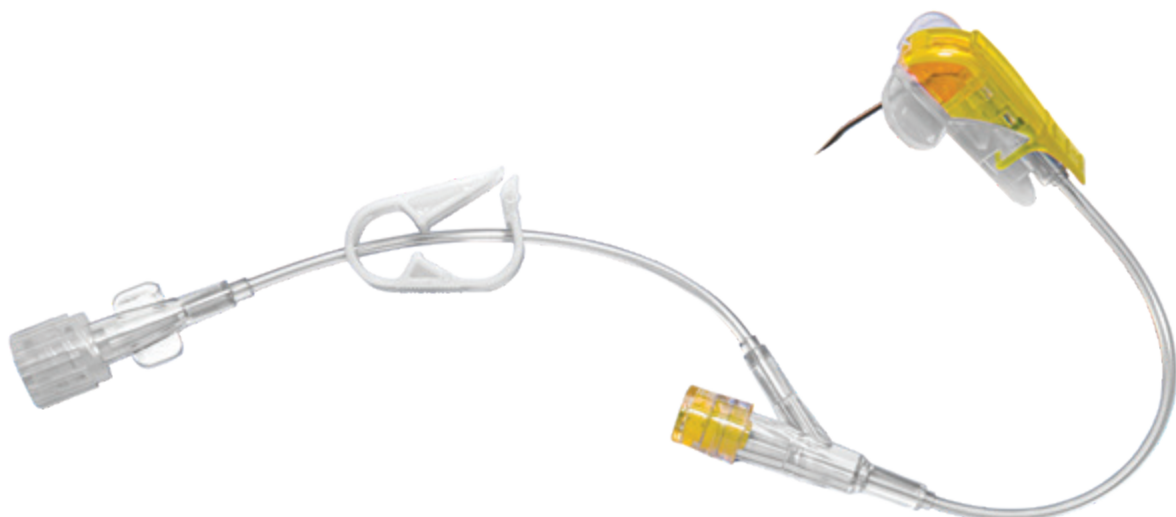

Safetouch Huber Needle set

STERILE, SINGLE-USE, NON-CORING HUBER NEEDLE INFUSION SET WITH SAFETY DEVICE

- Intended for the infusion of drug solutions into or blood sampling from a reservoir implanted in the body
- Each gauge and length of needle referenced above is power-injectable up to 300 psi
- Compact design without the need for additional stabilizing components
- Intuitive activation of the safety mechanism
- Designed for use with a vascular access infusion system
- Non-coring needle with cap
- Needle fixed at 90° angle
- Fixed colored wings
- Not made with latex or DEHP
- Available with or without Y-site
- Tubing lengths of 200 mm (with Y-site) and 200 or 450 mm (without Y-site)
- Available in 19, 20 and 22G
- Available in 19, 25, and 38 mm (3/4, 1, 1 1/2 inches respectively)

CE 0123

class IIb Medical Devices

Rule 7

MDD 93/ 42/ EEC

UMDNS 12729

Product Compliance

- Complies with the following norms, directives, and regulations:
 - ISO 9626:2016
 - ISO 7864
 - ISO 8536-4
 - ISO 10555-3
 - ISO 80369-7

Manufacturing DETAILS

Legal manufacturer: Nipro Corporation

Country of origin: Thailand

Sterilization and Shelf Life

Sterilization

EtO (Ethylene oxide)



Single-use only



Indicator sticker

Each outer box (or shipping carton) contains a chemical indicator that indicates sterility. This indicator comes in the form of a violet sticker that turns to blue when sterilized.


Shelf life

3 years



Materials Used

Not made with
 Latex or DEHP

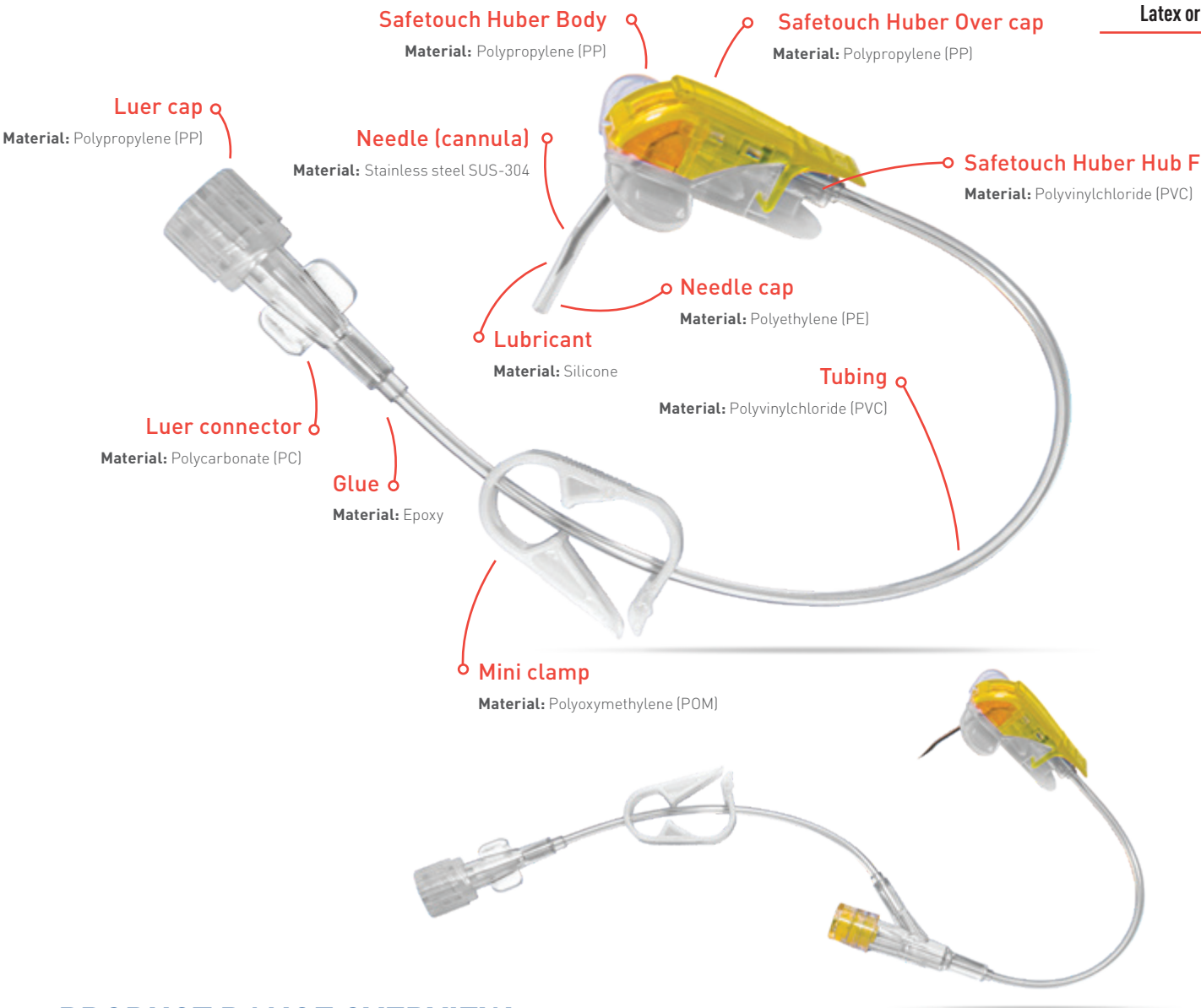

PRODUCT RANGE OVERVIEW

Reference code	Description	Needle outer diameter (G)	Needle length (mm)	Tube length volume (mm)	Y-site	Priming volume (ml)
HCNS-SF-1919-Y-ET	SF Huber coreless Ndl Set 19G 19 mm Y	19	19	200	Yes	0,6
HCNS-SF-1925-Y-ET	SF Huber coreless Ndl Set 19G 25 mm Y	19	25	200	Yes	0,6
HCNS-SF-2019-Y-ET	SF Huber coreless Ndl Set 20G 19 mm Y	20	19	200	Yes	0,6
HCNS-SF-2025-Y-ET	SF Huber coreless Ndl Set 20G 25 mm Y	20	25	200	Yes	0,6
HCNS-SF-2038-Y-ET	SF Huber coreless Ndl Set 20G 38 mm Y	20	38	200	Yes	0,6
HCNS-SF-2219-Y-ET	SF Huber coreless Ndl Set 22G 19 mm Y	22	19	200	Yes	0,6
HCNS-SF-2225-Y-ET	SF Huber coreless Ndl Set 22G 25 mm Y	22	25	200	Yes	0,6
HCNS-SF-2238-Y-ET	SF Huber coreless Ndl Set 22G 38 mm Y	22	38	200	Yes	0,6
HCNS-SF-1919-20-ET	SF Huber coreless Ndl Set 19G 19 mm 20 cm	19	19	200	No	0,3
HCNS-SF-1925-20-ET	SF Huber coreless Ndl Set 19G 25 mm 20 cm	19	25	200	No	0,3
HCNS-SF-1938-20-ET	SF Huber Coreless Ndl Set 19G 38 mm 20 cm	19	38	200	No	0,3
HCNS-SF-2019-20-ET	SF Huber coreless Ndl Set 20G 19 mm 20 cm	20	19	200	No	0,3
HCNS-SF-2025-20-ET	SF Huber coreless Ndl Set 20G 25 mm 20 cm	20	25	200	No	0,3
HCNS-SF-2038-20-ET	SF Huber coreless Ndl Set 20G 38 mm 20 cm	20	38	200	No	0,3
HCNS-SF-2219-20-ET	SF Huber coreless Ndl Set 22G 19 mm 20 cm	22	19	200	No	0,3
HCNS-SF-2225-20-ET	SF Huber coreless Ndl Set 22G 25 mm 20 cm	22	25	200	No	0,3
HCNS-SF-2238-20-ET	SF Huber coreless Ndl Set 22G 38 mm 20 cm	22	38	200	No	0,3

*Can be produced also in 450 mm on demand

Packaging Details

Without Y-site

Individual unit			Inner box				Outer box					
L (mm)	W (mm)	Weight (g)	Quantity (pcs)	L (mm)	W (mm)	H (mm)	Weight (kg)	Quantity (pcs)	L (mm)	W (mm)	H (mm)	Weight (kg)
98	85	13,2	20	308	105	152	0,376	160	434	322	332	3,77

With Y-site

Individual unit			Inner box				Outer box					
L (mm)	W (mm)	Weight (g)	Quantity (pcs)	L (mm)	W (mm)	H (mm)	Weight (kg)	Quantity (pcs)	L (mm)	W (mm)	H (mm)	Weight (kg)
98	85	13,8	20	308	105	152	0,388	160	434	322	332	3,868

Languages

26 languages: English (EN), Bulgarian (BG), Czech (CS), Danish (DA), German (DE), Greek (EL), Spanish (ES), Estonian (ET), Finnish (FI), French (FR), Croatian (HR), Hungarian (HU), Italian (IT), Lithuanian (LT), Latvian (LV), Dutch (NL), Norwegian (NO), Polish (PL), Portuguese (PT), Romanian (RO), Russian (RU), Slovakian (SK), Slovenian (SL), Serbian (SR), Swedish (SV), Turkish (TR).

Transport and storage conditions

Keep dry and away from sunlight.

Label details

	Individual unit	Inner box	Outer box
Product name	✓	✓	✓
Product size	✓	✓	✓
Quantity (pcs)	—	✓	✓
REF	✓	—	—
LOT	✓	✓	✓
	✓	✓	✓
STERILE EO	✓	✓	✓
CE	✓	✓	✓
	✓	✓	✓
	✓	✓	✓
	✓	✓	✓
	✓	✓	✓
EC REP	✓	✓	✓

	Individual unit	Inner box	Outer box
	—	✓	✓
	—	—	✓
	—	—	✓
	—	—	✓
	✓	✓	✓
Priming vol.	✓	—	—
MADE IN THAILAND	✓	✓	✓
Not made with LATEX	✓	✓	✓
Not made with DEHP	✓	✓	✓
	—	—	✓
Label/sticker area	—	✓	✓

