




AQUALiner[®]

Hydrophilic Nitinol Guidewire

Innovative coating with kink-resistant core and optimal tip design

Nipro's second-generation AQUALiner® Nitinol guidewire, with innovative hydrophilic coating and soft tip tapered design, is engineered to provide optimal pushability, trackability, and kink resistance in the most challenging anatomy while simultaneously providing lasting lubricity and tactile feel – helping you achieve procedural success.

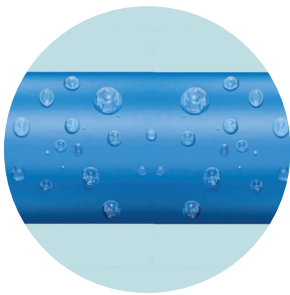


Innovative coating – Advanced coating technology produces a layer of moisture between the wire surface and either vessel walls or instrumentation

- Promotes smooth, frictionless navigation and maintains optimal tactile feel
- Maintains lubricity through multiple passes with minimal swelling



Kink-resistant navigation, super-elastic technology – The second-generation Nitinol (Ni-Ti) alloy core allows for the right blend of suppleness and rigidity that absorbs strain while maintaining durability

- Smooth, kink-free navigation through even the most tortuous vasculature
- 1:1 torque ratio for enhanced wire control and steerability

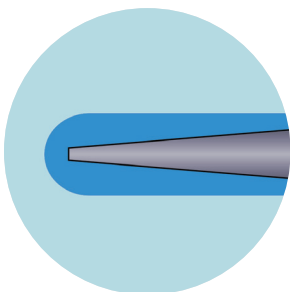

Soft tip design – A transitioned tip provides optimal stability while maintaining a flexible soft taper

- Soft tip for smooth, atraumatic insertion and navigation through vessels
- Minimal force is required to bend the tip

Control and durability

Angled tip

Straight tip

Available in stiff or regular wire
in various diameters and lengths,
packaged 5/box



Make AQUALiner® your go-to choice for smooth and easy navigation through tortuous vasculature.

Ordering information

STANDARD ORDER					
DIAMETER	LENGTH	TIP & WIRE COMBINATIONS			
		Straight-Regular	Straight-Stiff	Angled-Regular	Angled-Stiff
0.035"	150 cm	ALS+35150SR	ALS+35150SS	ALS+35150AR	ALS+35150AS
	180 cm	ALS+35180SR	ALS+35180SS	ALS+35180AR	ALS+35180AS
	260 cm	ALS+35260SR	ALS+35260SS	ALS+35260AR	ALS+35260AS
SPECIAL ORDER					
DIAMETER	LENGTH	TIP & WIRE COMBINATIONS			
		Straight-Regular	Straight-Stiff	Angled-Regular	Angled-Stiff
0.032"	150 cm	ALS+32150SR	-	ALS+32150AR	-
	180 cm	-	-	-	-
	260 cm	ALS+32260SR	-	ALS+32260AR	-
0.038"	150 cm	ALS+38150SR	ALS+38150SS	ALS+38150AR	ALS+38150AS
	180 cm	ALS+38180SR	ALS+38180SS	ALS+38180AR	ALS+38180AS
	260 cm	ALS+38260SR	ALS+38260SS	ALS+38260AR	ALS+38260AS

Straight-Regular and Angled-Regular = blue label color. Straight-Stiff and Angled-Stiff = brown label color.

AQUALiner[®]

Hydrophilic Nitinol Guidewire

YOUR LINK BETWEEN PATIENT AND CARE

Nipro Medical Corporation is a subsidiary of Nipro Corporation Japan, established in 1954 and a world leader in healthcare with over 35,000 employees around the globe. We are committed to delivering innovative, high quality solutions in the form of renal, medical-surgical, and vascular products. This fundamental objective drives us to develop cutting-edge, value-added products for healthcare professionals that are safe, effective, and easy to use.

As an integral part of a vast global manufacturing network, Nipro Medical Corporation's North America division is responsible for managing the sales, marketing, and business operations of the United States and Canada – unified in our commitment to improving patient care and quality of life.



NIPRO MEDICAL CORPORATION

200 Crossing Blvd., Bridgewater, NJ 08807
T: +1 (908) 393-7030
www.nipro.com | info@nipromed.com

© 2022 Nipro Medical Corp. All Rights Reserved.
Nipro and AQUALiner are registered trademarks of Nipro Corporation. Printed in U.S.A. RA-BRUSAQUA Rev. 3
NMC-LBL-0007_05 May 2022-USA-EN