



CAPTUS[®]

Vascular Retrieval System



www.nipro.com

CAPTUS®

Vascular Retrieval System

For simple & complex retrievals, reach for one device

CAPTUS provides the optimal balance of stiffness, flexibility, and precision in one easy to use vascular retrieval device.

Thanks to two pairs of independent suture loops and a unique nitinol braided frame, CAPTUS offers a secure and confident capture, enabling you to be prepared for filter retrieval complications along the way.



1. Self-centering frame design

A unique braided frame design self-centers the hook and/or apex of the filter into the loop snare, **allowing for a straightforward and secure hook capture.**

2. Braided loop snare


A proprietary crosshair, inner loop suture design enables a nearly infinite number of snare points, **simplifying capture and minimizing filter release.**

3. Shapeable PEEK shaft

A shapeable polyether-ether-ketone (PEEK) shaft, **allowing you to shape the CAPTUS shaft** to direct the self-centering cone toward filter hooks that may not be centered in the inferior vena cava (IVC).

4. Physician-tensioning hub

A precise, physician-tensioning feature tightens the loop snare around the filter hook and/or apex, **ensuring a secure grasp for a confident capture.**

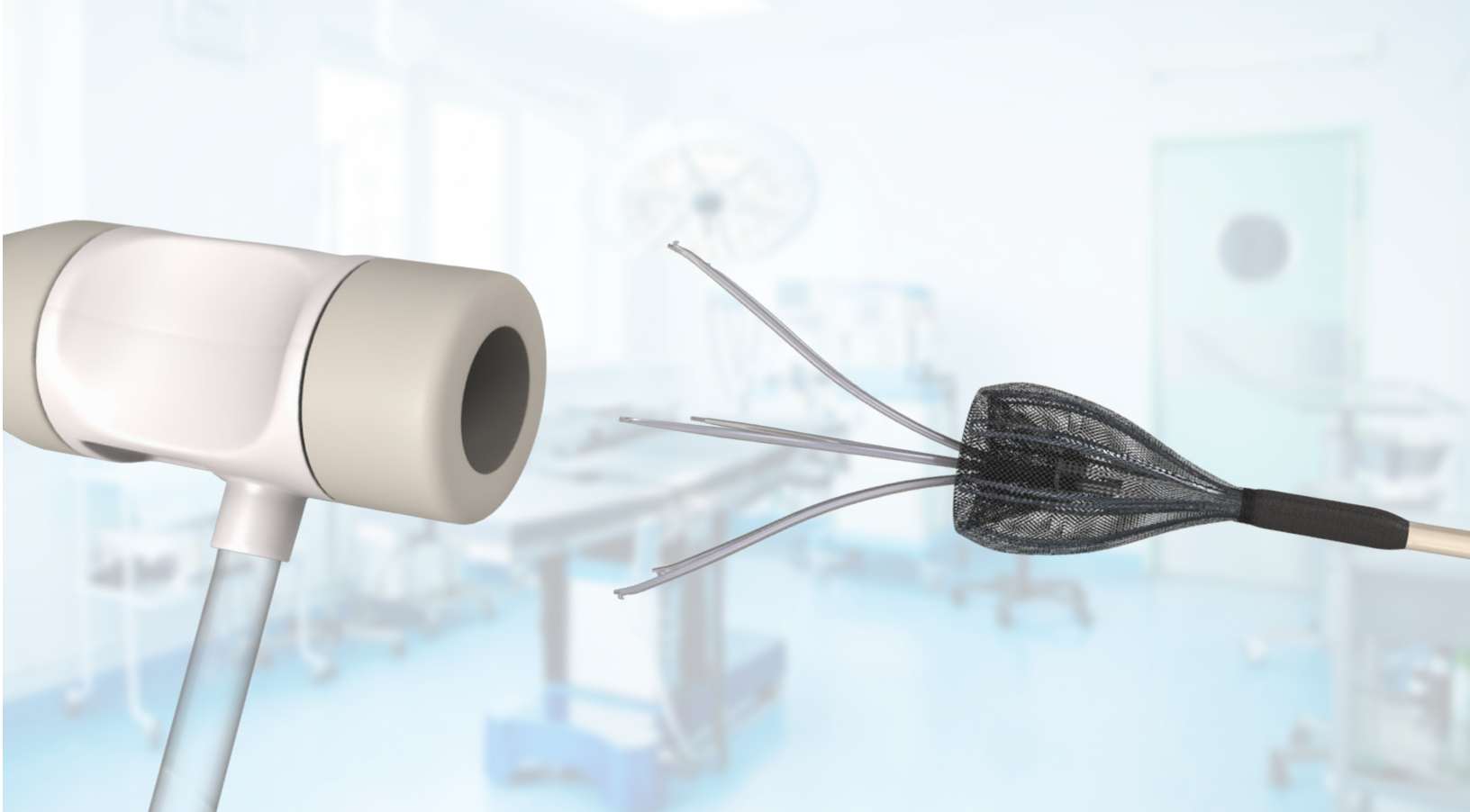


DON'T GET CAUGHT OFF GUARD

Make **CAPTUS** your go-to choice, preparing you for routine retrievals that unexpectedly go from simple to complicated.

Specifications

Order number	Item	Size
CAP1412	CAPTUS Vascular Retrieval Device	14 mm braided loop snare on 120 cm PEEK shaft with physician-tensioning feature

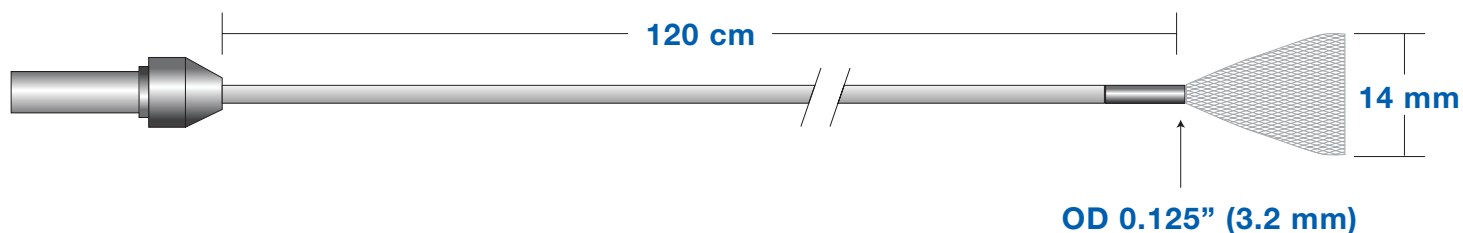

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. The CAPTUS Vascular Retrieval System is intended for use in the cardiovascular system to retrieve foreign objects, including, but not limited to guidewires, coils, balloons, catheters, and filters.

CAPTUS[®]

Vascular Retrieval System

YOUR LINK BETWEEN PATIENT AND CARE

Nipro Medical Corporation is a subsidiary of Nipro Corporation Japan, established in 1954 and a world-leader in healthcare with over 30,000 employees around the globe. We are committed to delivering innovative, high quality solutions in the form of renal, medical-surgical, and vascular products. This fundamental objective drives us to develop cutting-edge, value-added products for healthcare professionals that are safe, effective, and easy to use.

As an integral part of a vast global manufacturing network, Nipro Medical Corporation's North America division is responsible for managing the sales, marketing, and business operations of the United States and Canada – unified in our commitment to improving patient care and quality of life.



NIPRO MEDICAL CORPORATION

200 Crossing Blvd., Bridgewater, NJ 08807
T: +1 (908) 393-7030
www.nipro.com | info@nipromed.com

© 2021 Nipro Medical Corp. All Rights Reserved.
CAPTUS[®] is a registered trademark of Avantec, A Nipro Company.
NMC-LBL-0060_13 October 2021-USA-EN